

Shriners Gait Model Setup with Multi-segment Foot Model (Currently Vicon Use Only)

1. Moving VST files to correct location for Vicon (Vicon Skeletal Models)

- a. Copy the .vst files from the vsts folder in the GitHub project (ViconClinicalGait/vsts/)
 - i. Current vsts included with the project are:
 1. SMACNet Foot Model no KAD.vst (Multi-segment Foot Model)
 2. SMACNet Foot Model no KAD FootNotFlat 2ndStep (Multi-segment Foot Model)
- b. Paste into your Vicon ModelTemplates folder
 - i. Example filepath
C:\Users\Public\Documents\Vicon\Nexus2.x\ModelTemplates


2. Moving Pipeline files to correct location for Vicon

- a. Copy the .Pipeline files from the Utils folder in the GitHub project (ViconClinicalGait/Utils/Pipeline)
 - i. Current Pipelines included with the project are:
 1. Py3_Static
 2. Py3_Dynamic
 3. Py3_SittingFoot_Static_BF
- b. Paste into your Vicon Pipelines folder
 - i. Example filepath
C:\Users\Public\Documents\Vicon\Nexus2.x\Configurations\Pipelines

3. Moving Python Code to the correct location for Vicon

- a. Copy the .py files from the Py3_ShrineGaitModel folder in the GitHub project (ViconClinicalGait/Py3_ShrineGaitModel)
 - i. Current Python files included with the project are:
 1. Py3_CreateGCD
 2. Py3_DynamicMain
 3. Py3_GaitModules
 4. Py3_MathModules
 5. Py3_StaticMain
 6. Py3_SittingFootStaticMain.py
 7. Py3_WriteHeelToeMarkers.py
 8. Py3_UserPreferences
- b. Paste into a Vicon Python folder
 - i. Example filepath: C:\Users\Public\Documents\Vicon\Python
- c. Possible errors in where these files are placed could impact their use
- d. File location matters significantly, the Py3_UserPreferences.py file must be in the same directory(folder) as the Py3_StaticMain.py file for it to identify the appropriate file.

4. **Updating the User Preferences File:** Several items need to be reviewed in “UserPreferences.py” for each individual lab, these are the contained in the first 7 lines of code
- a. The marker diameter needs to be entered for your lab
 - i. Most likely the diameter will be 9.5 mm
 - ii. If you don’t use a base that adds thickness the values entered will be:
 1. `self.MarkerDiameter = 9.5`
 - iii. If you use a base, add the thickness of the base and multiply by 2. For example if you use a 2 mm base add the value of 4 as follows:
 1. `self.MarkerDiameter = 9.5+4`
 - b. For the next line of code (`self.StaticForwardDirection =`) you will need to know the global axis coordinate system for your lab.
 - i. Determine the global axis coordinate system
 - ii. Determine the direction that participants stand during the static trial (+X, -X, +Y, -Y). Enter the direction the participant stands with respect to the lab coordinate system:
 - iii. Here are 4 examples:

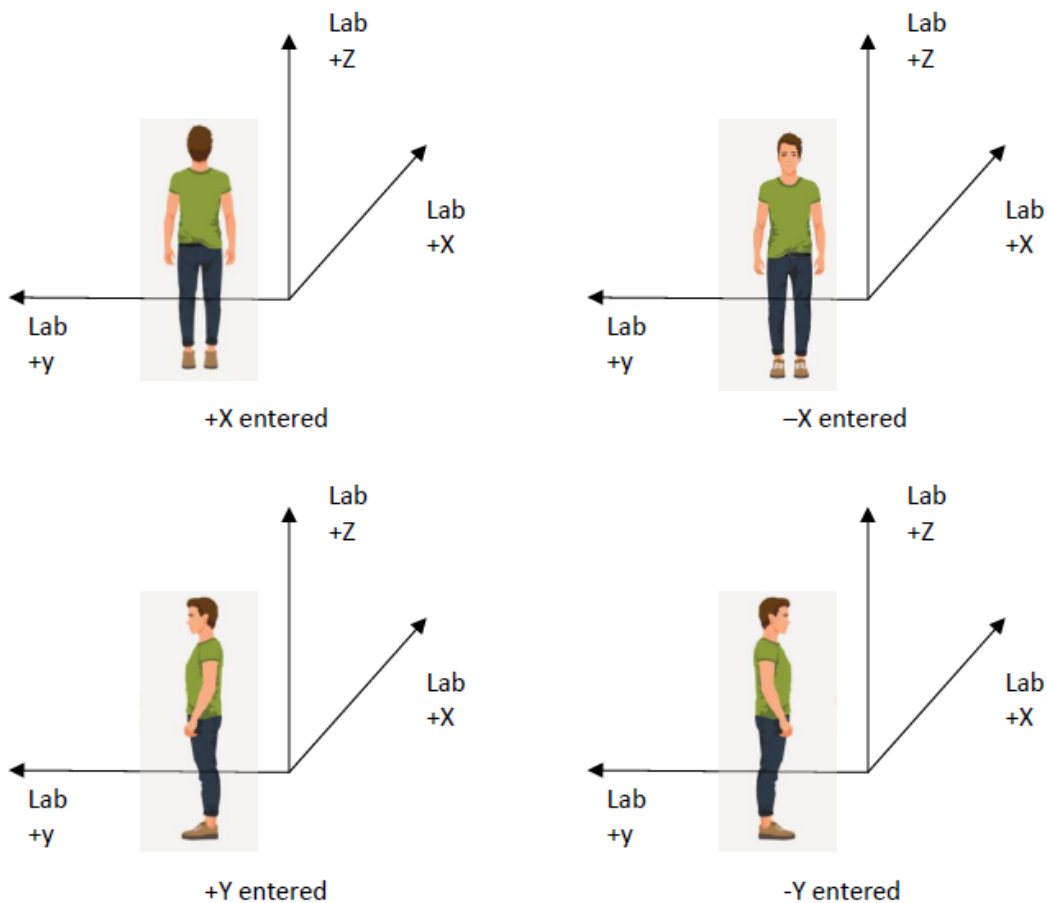

- iv. Enter the correct '+X' or '-X' or '+Y' or '-Y' for your lab and data collection protocol: e.g. `self.StaticForwardDirection = '-X'`
- c. The 3rd line of code is related to the selection of the hip model implementation (`self.HipModelName`)
 - i. The default for clinical gait analysis for Shriners Children's is 'HarringtonHybrid'
 - ii. 'HarringtonHybrid' uses regression equations to determine the x,y,z coordinates of the hip joint centers. The full set of equations use Pelvic Width (measured distance between ASIS during physical exam), Pelvic Depth (distance between Midpoint of ASIS markers and Midpoint of PSIS markers), and Leg Length from physical exam. BOTTOM LINE: leg length and ASIS-width need to be measured on physical exam and used.
- d. The 4th and 5th lines establish the rotation sequence for the pelvis and trunk. The 2 options are the orders ROT (Rotation-Obliquity-Tilt) or TOR (Tilt-Obliquity-Rotation). **'ROT' is the setting that should be used.**
 - i. `self.TrunkRotationSequence = 'ROT'`
 - ii. `self.PelvisRotationSequence = 'ROT'`
- e. The 6th line sets the shank coordinate system as Distal or Proximal. **The setting of 'Distal' should be used.** NOTE: Distal or Proximal must start with capital letter
 - i. `self.ShankCoordinateSystem = Distal`
- f. Lines 7 to 70 establish the marker labels used in the VST files. The labels in the `py3_UserPreferences.py` file match these labels as needed, so this part of the file does not need to be updated unless you decide to change the marker labels.

5. Static Trial

- a. In session folder, create subject and attach the appropriate vst (SMACNet Foot Model no KAD)
 - i. Under Subject Properties, enter all required anthropometrics (outlined in red)
- b. Record a static trial
 - i. Open static
 - ii. Run "Reconstruct and Label" pipeline. Check marker labels

General	
Bodymass (kg):	
Height (mm):	
InterASISDistance (mm):	
Left	
LegLength (mm):	
ASISTrochanterDistance (mm):	
KneeWidth (mm):	
AnkleWidth (mm):	
VarValAngle (deg):	
CalcanealPitch (deg):	
HindfootProgression_relBiMal (deg):	
1stRayPitch (deg):	
ForefootProgression_relBiMal (deg):	
Right	
LegLength (mm):	
ASISTrochanterDistance (mm):	
KneeWidth (mm):	
AnkleWidth (mm):	
VarValAngle (deg):	
CalcanealPitch (deg):	
HindfootProgression_relBiMal (deg):	
1stRayPitch (deg):	
ForefootProgression_relBiMal (deg):	

c. Run “Py_3 Static” pipeline

i. On the following screen:

1. Fill out required patient/subject information (outlined in red)
2. Match Static Frame Number to the selected frame number of the static trial (outlined in blue). Static frame should be fully labeled.
3. Check that all required anthropometric measures are filled in (outlined in green)
4. Foot Model Specific Parameters
 - a. Under “Foot Model” (outlined in purple), check the boxes for both feet
 - b. Make sure Hindfoot Varus/Valgus is populated
 - c. Other measures are optional (we highly recommend reviewing the **4.FootModelSOP** if using the foot model to help you determine if you would like to use additional xray measures)

Static

Patient Information

First Name Last Name

Patient Number

Date Of Birth 01 Jan 1981

Trial Activity Static

Trial Modifier Barefoot

Static Forward Direction +X

Assistive Device None

Anthropometric Parameters

Static Frame Number 85

Static File Anthropometric File Static_BF_2025Consistency

Date Of Data Capture 17 Dec 2025

kg 80.0 BodyMass

mm 1800 Height

mm 248 ASIS-to-ASIS Distance

mm 995 995 Leg Length

mm 127 126 Knee Width

mm 78 78 Ankle Width

mm 80 80 ASIS-to-GT Distance

Foot Model

☒ Left ☒ Right

deg [-Val/+Var] 0 Hindfoot Varus/Valgus

deg [+Up/-Down] Calcaneal Pitch

deg [+Int/-Ext] Hindfoot Progression [rel. to Bimalleolar axis]

deg [-Up/+Down] First Metatarsal Pitch

deg [+Int/-Ext] Forefoot Progression [rel. to Bimalleolar axis]

Special Test Conditions

☒ Left ☒ Right Subject Plantigrade during Static Trial

☐ Subject Shod (with or without orthoses)

mm Left Right Difference in sole thickness between heel and toes (mm)

☒ N/A ☐ Iliac ☐ Triad Apply Pelvic Fix Option

Knee Option ☐ KAD ☒ M/L Markers

Save & Proceed

Quit

5. Under “Special Test Conditions”, select appropriate boxes to indicate whether subject is plantigrade.
 - a. Note: If using the foot model as is, the subject **MUST** be plantigrade during static collection
 - i. If they cannot stand plantigrade, please see **4. Foot Model SOP** (Appendix A) for information on how to perform a sitting static where the feet are plantigrade
 - ii. The appendix follows a similar flow to this user guide, with detailed instructions on how to collect each trial and which options to select.
6. We recommend not changing the knee option setting. This model will work with a KAD (knee alignment device), but we have shifted to using medial and lateral knee markers
7. Ignore the Pelfix option (a separate user guide is included for using this option)
8. Click “Save & Proceed”
- ii. Compare marker measurements to clinical measurements (The Shriners Gait Model uses clinical measurements to calculate joint centers and virtual markers). Marker measurements should be relatively close to clinical measures in most cases (5-10 mm for knees and ankles, around 20-25 mm for hip width). More significant differences should be evaluated to determine if marker placement is in line with the clinical measurement. Click “Proceed”
- iii. Check Standing Posture measurements. If you are using the Standard Model, the Foot Segment measurements will be empty. Click “Save Results” and “Save PDF”. Click “Quit”
 1. “Save Results” saves a static.py file that contains all of the static trial information needed to process the dynamic trials (marker relationships between technical and anatomical coordinate systems, anthropometrics, etc.)
 2. “Save PDF” saves a pdf output report of the static trial you just processed, showing the marker measurement comparison, standing posture measurements, etc.


6. Dynamic Trials

- a. Medial knee and medial malleoli markers can be removed for dynamic trials
- b. Record dynamic walking trials in your capture volume

7. Processing Dynamic Trials

- a. Open dynamic walking trial
- b. Run “Reconstruct & Label” pipeline
 - i. Check marker labels and gap fill as needed

- ii. Crop trial as needed. The Shriners Gait Model requires all markers to be present on the first frame of the trial
- c. Add in gait events
 - i. Label at least one full gait cycle as shown below

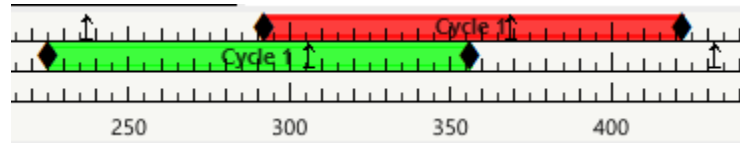

- d. This model does support kinetics requiring force plates to be integrated with Vicon
 - i. Assign each force plate based on which foot strikes the plate (right or left)
- e. Run "Py3_Dynamic" pipeline

8. Outputs

- a. The Py3_Dynamic pipeline will complete the calculations and save the results (kinematics, kinetics, etc.) as model outputs
- b. These outputs can be found in Vicon by:
 - i. clicking on the Subjects Tab
 - ii. clicking the drop down menu to the left of the subject name
 - iii. clicking the drop down menu to the left of Model Outputs
 - iv. and clicking the drop down menu to the left of the model output you are interested in
- c. The pipeline will also output a .GCD file
 - i. This is a type of text file that contains each model output time normalized to be represented as percentage of the gait cycle (101 points 0-100%)
 - ii. This text file can be used as a coding input to visualize the data
 - iii. No visualization code has been provided